

GC Vitamin D Binding Protein Antibody

Mouse Monoclonal Antibody [Clone VDBP/4481]


Catalog No	Format	Size
2638-MSM1-P0	Purified Ab with BSA and Azide at 200ug/ml	20 ug
2638-MSM1-P1	Purified Ab with BSA and Azide at 200ug/ml	100 ug
2638-MSM1-P1ABX	Purified Ab WITHOUT BSA and Azide at 1.0mg/ml	100 ug
Applications	Tested Dillution Note	

Applications	Tested Dillution	Note
Immunohistochemistry (IHC)	5	30 min at RT. Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95°C followed by cooling at RT for 20 minutes


Product Details	
Clone	VDBP/4481
Gene Name	GC
Immunogen	Recombinant fragment (around aa 35-175) of human GC protein (exact sequence is proprietary)
Host	Mouse
Clonality	Monoclonal
Isotype / Light Chain	IgG2b / Kappa
Mol. Weight of Antigen	53kDa.
Cellular Localization	Secreted
Species Reactivity	Human
Positive Control	Human liver, kidney or pancreas.

*Optimal dilution for a specific application should be determined.

Product Images for GC Vitamin D Binding Protein Antibody

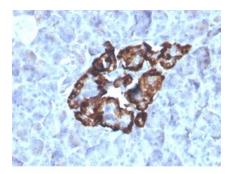


IHC analysis of formalin-fixed, paraffin-embedded human renal oncocytoma. Staining using VDBP/4481 at 2ug/ml in PBS for 30min RT. HIER: Tris/EDTA, pH9.0, 45min. 2°C: HRP-polymer, 30min. DAB, 5min.

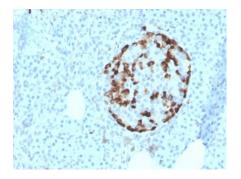

IHC analysis of formalin-fixed, paraffin-embedded human pancreas. Staining using VDBP/4481 at 2ug/ml in PBS for 30min RT.




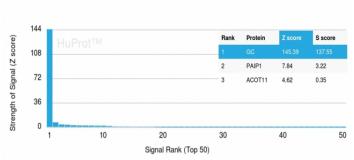

Analysis of Protein Array containing more than 19,000 full-length human proteinsusing GC Vitamin D Binding Protein Mouse Monoclonal Antibody (VDBP/4481). Z- and S-Score: The Z-score represents the strength of a signal that a monoclonal antibody (MAb) (in combination with a fluorescently-tagged anti-IgG secondary antibody) produces when binding to a particular protein on the HuProTM array. Z-scores are described in units of standard deviations (SD's) above the mean value of all signals generated on that array. If targets on HuProtTM are arranged in descending order of the Z-score, the S-score is the difference (also in units of SD's) between the Z-score. Sscore therefore represents the relative target specificity of a MAb to its intended target. A MAb is considered to specific to its intended target, if the MAb has an S-score of at least 2.5. For example, if a MAb binds to protein X with a Z-score of 43 and to protein Y with a Z-score of 14, then the S-score for the binding of that MAb to protein X is equal to 29.

Specificity & Comments

Vitamin D-binding protein (DBP) is a multi-functional serum protein that binds to the plasma membranes of numerous cell types and mediates a variety of cellular functions. The locus of the DBP protein (also known as group-specific component protein or GC) is located at human chromosome 4q13.3. DBP functions in organ-specific transportation of vitamin D and its metabolites to the various target organs of the vitamin D endocrine system. In addition, DBP has immunomodulatory properties and is able to bind to the surface of leukocytes. DBP binds to the plasma membrane through a chondroitin sulfate proteoglycan. DBP serves as a co-chemotactic factor for C5a to enhance the chemotactic activity of C5a. DBP can also bind to globular Actin with high affinity and is involved in the clearance of Actin from the blood. DBP plays an important role in osteoclast differentiation. The diverse cellular functions of DBP require its cell surface binding ability to mediate different biological processes.

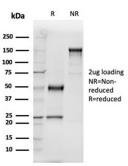

SDS-PAGE Analysis of Purified GC Mouse Monoclonal Antibody (VDBP/4481). Confirmation of Purity and Integrity of Antibody.

Supplied As

200ug/ml of Ab purified from Bioreactor Concentrate by Protein A/G. Prepared in 10mM PBS with 0.05% BSA & 0.05% azide. Also available WITHOUT BSA & azide at 1.0mg/ml.

Storage and Stability

Antibody with azide - store at 2 to 8 °C. Antibody without azide - store at -20 to -80 °C. Antibody is stable for 24 months. Non-hazardous. No MSDS required.

Research Areas

Cardiovascular

Limitations and Warranty

This antibody is available for research use only and is not approved for use in diagnosis. There are no warranties, expressed or implied, which extend beyond this description. Company is not liable for any personal injury or economic loss resulting from this product.

