

CD47 / IAP (Integrin Associated Protein) Antibody

Mouse Monoclonal Antibody [Clone CD47/3019]

Catalog No	Format	Size
961-MSM9-P0	Purified Ab with BSA and Azide at 200ug/ml	20 ug
961-MSM9-P1	Purified Ab with BSA and Azide at 200ug/ml	100 ug
961-MSM9-P1ABX	Purified Ab WITHOUT BSA and Azide at 1.0mg/ml	100 ug


Applications	Tested Dillution	Note
Flow Cytometry (Flow)	1-2ug/million cells	
Immunofluorescence (IF)	1-3ug/ml	
Immunohistochemistry (IHC)	1-2ug/ml	30 min at RT. Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95°C followed by cooling at RT for 20 minutes

Product Details

Clone	CD47/3019
Gene Name	CD47
Immunogen	A recombinant fragment (around aa 18-135) of human CD47 protein (exact sequence is proprietary)
Host	Mouse
Clonality	Monoclonal
Isotype / Light Chain	IgG1 / Kappa
Mol. Weight of Antigen	~50kDa
Cellular Localization	Cell membrane, Cell surface
Species Reactivity	Human
Positive Control	MCF-7
*Optimal dilution for a specific app	plication should be determined.

dilution for a specific applic ation snould be determin

Product Images for CD47 / IAP (Integrin Associated Protein) Antibody

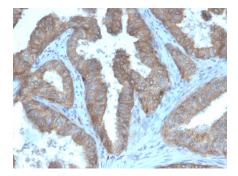


Formalin-fixed, paraffin-embedded human prostate carcinoma stained with CD47 Mouse Monoclonal Antibody (IAP/3019).


kDa		R	NR	
250				
150			-	
100				
75	-			2ug loading
50	-	-		NR=Non- reduced
37	-			R=reduced
25	_	_		
20				
15	_			
10				

SDS-PAGE Analysis Purified CD47 Mouse Monoclonal Antibody (IAP/3019). Confirmation of Integrity and Purity of Antibody.



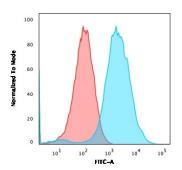


Flow Cytometric Analysis of PFA-fixed MCF-7 cells using CD47 Mouse Monoclonal Antibody (IAP/3019) followed by goat anti-mouse IgG-CF488 (blue); isotype control (red).

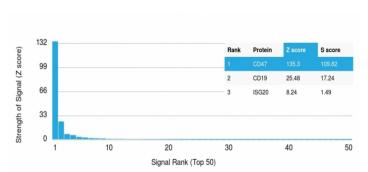

Analysis of Protein Array containing more than 19,000 full-length human proteinsusing CD47 Mouse Monoclonal Antibody (CD47/3019) Z- and S- Score: The Z-score represents the strength of a signal that a monoclonal antibody (MAb) (in combination with a fluorescently-tagged anti-IgG secondary antibody) produces when binding to a particular protein on the HuProtTM array. Z-scores are described in units of standard deviations (SD's) above the mean value of all signals generated on that array. If targets on HuProtTM are arranged in descending order of the Z-score, the S-score is the difference (also in units of SD's) between the Z-score. S-score therefore represents the relative target specificity of a MAb to its intended target. A MAb is considered to specific to its intended target, if the MAb has an S-score of at least 2.5. For example, if a MAb binds to protein X with a Z-score of 43 and to protein X with a Z-score of 14, then the S-score for the binding of that MAb to protein X is equal to 29.

Immunofluorescence staining of PFA-fixed MCF-7 cells. CD47 Mouse Monoclonal Antibody (IAP/3019) followed by goat anti-mouse IgG-CF488 (green), Nuclei are stained with RedDot.

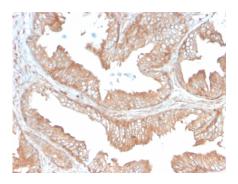

Formalin-fixed, paraffin-embedded human prostate carcinoma stained with CD47 Mouse Monoclonal Antibody (IAP/3019).

Specificity & Comments

This antibody reacts with Ig domain of CD47 protein. CD47, originally named integrin-associated protein (IAP), is a 50kDa protein containing five membrane-spanning sequences and a short cytoplasmic tail. CD47 plays a role in both cell adhesion by acting as an adhesion receptor for THBS1 on platelets, and in the modulation of integrins. It is important in memory formation and synaptic plasticity in the hippocampus. CD47 may play a role in membrane transport and/or integrin dependent signal transduction.

Supplied As

200ug/ml of Ab purified from Bioreactor Concentrate by Protein A/G. Prepared in 10mM PBS with 0.05% BSA & 0.05% azide. Also available WITHOUT BSA & azide at 1.0mg/ml.

Storage and Stability

Antibody with azide - store at 2 to 8 °C. Antibody without azide - store at -20 to -80 °C. Antibody is stable for 24 months. Non-hazardous. No MSDS required.

Research Areas

Cardiovascular, Immunology, Immune checkpoint

Limitations and Warranty

This antibody is available for research use only and is not approved for use in diagnosis. There are no warranties, expressed or implied, which extend beyond this description. Company is not liable for any personal injury or economic loss resulting from this product.

